

## Seers Treatment Tables

Seers Medical, an innovative medical engineering company with over 20 years of experience manufacturing Akron tables, is the leading UK manufacturer of physiotherapy, rehabilitation and medical tables. Now offering a broad line of treatment tables, in the following pages we are showing you the most popular models amongst Canadians physiotherapists. These models have been engineered with the clinicians' needs in mind. Seers distributes to over 60 countries worldwide; OrthoCanada is very fortunate and proud to be the exclusive Canadian distributor for the Seers Medical product line.



### Anti-Entrapment Gas Struts

The backrest and head section are fitted with power-assisted anti-entrapment gas struts which provide a very safe and reliable control mechanism. They require the use of the release lever only when returning the section back to the "closed or negative" position. This unique design greatly reduces the possibility of the gas strut seals being damaged by misuse or with accidents.

### Self Lubricating Bearings

Precision made self lubricating bearings which ensure that the couch frame repeatedly raises and lowers in a smooth and silent manner with minimal maintenance.

### Sturdy Wheels

The wheels are 3" in diameter and have a combined load capacity of 660lbs which is more than enough to support the weight of the table with a patient and therefore suitable for many years of use. Not to be used to transport patients.

### Anti-Bacterial Upholstery

1.6" medium density foam, covered with a waterproof, anti-bacterial, fire retardant Mundial vinyl. The anti-bacterial impregnated upholstery protects against microbial and fungal growth including salmonella and E coli. This vinyl has a very high Martindale rating of 100,000 rubs which indicates that it's extremely hard wearing and will last for many years if treated and cleaned correctly. (Martindale is an abrasion test and most vinyls are rated at 50,000).

### One Hand Release

One hand release head section. Fitted with power assisted anti-entrapment gas struts. (Refer to description on opposite page).

### Treated Frames with High Load Capacity

All table frames are treated against rust prior to being coated, and stove enamelled with epoxy powder paint which includes an anti-microbial additive for additional infection control and ease of cleaning. 3mm thick steel walls to provide extreme strength and product longevity. Safe lifting capacity of 528lbs and the capacity to be upgraded to 715lbs purely by fitting a more powerful motor. Tables are tested with 2.5 times their lifting capacity (tested with 1320 lbs).

### Retractable Wheel Design

Activated at each end of the table. Simply push the lever down to raise the table on its wheels for easy transport. Once the table is properly positioned, flick the lever upward to drop the table down onto its 4 legs. This wheel system is a well proven design and has been produced for over 30 years.


### Perimer Foot Switch

This is the best design, by far. The Seers perimeter foot switch is motionless. No need to push it up/down or from side to side to activate the motor. Simple foot pressure on the top or bottom of the switch membrane will activate the table height upwards or downwards as required, leaving the clinicians' hands free to support or work with their patient. This design has proven to be trouble free.

### Adjustable Leveling Foot

Adjustable leveling foot for uneven flooring can easily be adjusted to set all 4 legs firmly on the ground for added stability.

### Stock Colors

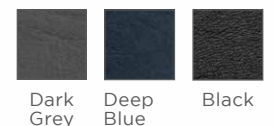

Dark Grey, Deep Blue, Black

### Custom colors (more available upon request)



Light Grey, Cream, Sky Blue, Lotus Green, Cherry, Burgundy

Mundial vinyl with marindale rating of 100 000 rubs