


# LYMPHASTIM APPLICATORS

## OVERLAPPING CHAMBERS TECHNOLOGY

All Lymphastim applicators feature overlapping chambers for precise and more efficient therapy.

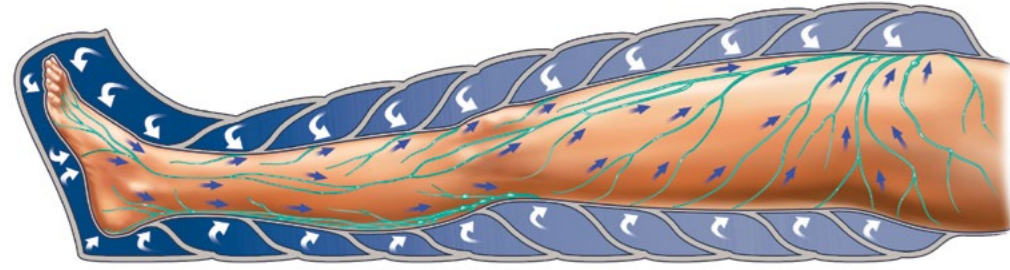

### FEATURES & BENEFITS OF THE BTL-6000 LYMPHASTIM APPLICATORS

- Ergonomically designed for maximum comfort.
- Connectors for easy applicator change by a single click.
- Automatic sleeve emptying after treatment for comfortable handling.
- Automatic identification of the applicators; the unit automatically identifies the connected applicator and select the compression flow accordingly.

### WIDE RANGE OF APPLICATORS

#### TROUSERS 24

The velcro applicator perfectly fits all body types due to the large coverage area of the velcro.



#### TROUSERS 24 SIZE XL

For extra large patients.



#### TROUSERS 24 WITH ZIPPER

Full-length, high-quality zipper for quick and easy dressing / undressing.



#### TROUSERS 24 UNISIZE

The unisize trousers have been designed to change size fast and easily.



#### LEG APPLICATORS

A selection of three sizes (S, M and L) is available.



#### ARM APPLICATORS

The arm applicators are unique due to their extremely ergonomical large breast chambers.



## BTL-6000

### LYMPHASTIM



sales@btlnet.com  
www.btlnet.com

All rights reserved. Although every care has been taken to provide accurate and up-to-date information, no responsibility can be accepted for any mistakes or misprints contained in this catalogue. The products and their specifications, as shown in this catalogue, may change without notice.

# BTL-6000 LYMPHASTIM UNRIVALLED SYSTEM FOR PRESSOTHERAPY

Lymphatic drainage is well-established & clinically-proven therapeutic method which brings immediate and visible results.

BTL's Lymphastim works on the pneumatic pressotherapy principle. Special applicator sleeves with multiple overlapping chambers provide gentle massage that encourages the natural circulation of the lymph through the body.

## DIFFERENT FIELDS OF APPLICATION

- Rehabilitation and sports medicine
- Lymphology and rheumatology clinics
- Spas, wellness hotels and beauty centres

## MOST COMMON INDICATIONS

### MEDICAL

- Chronic oedemas
- Post-surgical obstruction of lymph vessels
- Venous insufficiency and ischemic diseases

### WELLNESS AND AESTHETICS

- Cellulite
- Postliposuction care
- Heavy-legs syndrome



## BTL-6000 LYMPHASTIM TOPLINE

### PROFESSIONAL 12-CHANNEL PRESSOTHERAPY DEVICE

- 15 programmes for various compression flow
- Built-in encyclopaedia with treatment recommendations
- Memory for up to 500 user-defined programmes
- Advanced setting for customized treatment

### FEATURES & BENEFITS OF THE BTL-6000 LYMPHASTIM TOPLINE

- 26 preset therapy protocols and sequences of programmes
- Individually adjustable pressure for each chamber
- Working chamber activity shown on display
- Option to disable one or more chambers where the patient cannot tolerate the pressure
- Advanced settings of pressure gradient

## BTL-6000 LYMPHASTIM EASY

### 12-CHANNEL / 6-CHANNEL ECONOMY MODELS WITH LED INTERFACE

- 8 preset programmes that can be modified by the user
- Basic pressure gradient and programme speed setting
- Pressure can be changed during the course of therapy
- Working chamber activity shown on bar graph



# TECHNICAL PARAMETERS

## TECHNICAL SPECIFICATIONS OF THE BTL-6000 LYMPHASTIM

MODEL	BTL-6000 LYMPHASTIM 12 TOPLINE	BTL-6000 LYMPHASTIM 12 EASY	BTL-6000 LYMPHASTIM 6 EASY
Ordering number	P6000.003	P6000.002	P6000.001
User-interface	5.7" touch screen	LED display	LED display
BTL Encyclopaedia	26 clinical protocols (programmes in sequences)		
User programmes	500 user-defined programmes		
Preset programmes	15	8	8
Numbers of channels	12	12	6
Pressure setting	20-160 mmHg	20-160 mmHg	20-160 mmHg
Gradient	0-100 % continuous setting	0, 10, 20, 30 %	0, 10, 20, 30 %
Dimensions	320 x 190 x 280 mm		
Weight	7,5 kg		
Mains supply	100-240 V, 50-60 Hz		
Standard accessories	Basic tube 12 Topline, power cord	Basic tube 12 Easy, power cord	Basic tube 6 Easy, power cord

### Optional accessories for Lymphastim 12 Easy / 12 Topline

Part number P6000.060	Trousers 24 chambers	Part number P6000.075	Arm 8 chambers-size S
Part number P6000.059	Trousers 24 chambers with velcro	Part number P6000.067	Leg 10 chambers-size L
Part number P6000.071	Trousers 24 chambers with velcro-size XL	Part number P6000.066	Leg 10 chambers-size M
Part number P6000.072	Trousers 24 chambers with zipper		
Part number P6000.065	Arm 8 chambers		

### Optional accessories for Lymphastim 6 Easy

Part number P6000.061	Trousers 12 chambers	Part number P6000.063	Leg 6 chambers-size L
Part number P6000.064	Arm 6 chambers	Part number P6000.062	Leg 6 chambers-size M
Part number P6000.074	Arm 6 chambers-size S	Part number P6000.077	Leg 6 chambers-size S

### Optional accessories for all units

Part number P6000.034	Interface for 2 arms or legs	Part number P6000.073	Extension strip for trousers 24 chambers with zipper
Part number P6000.070	Extension strip for trousers 24, 12	Part number BTL-13.026	BTL-1300 Lymphatic couch
Part number P6000.069	Extension strip for leg applicators-size L	Part number T013.025	BTL-1300 Lymphatic couch Basic
Part number P6000.068	Extension strip for leg applicators-size S,M		
Part number P5002.009	Trolley		
Part number P6000.211	Plastic trolley		

